


CircAid[®]

MEDICAL PRODUCTS

The Graduate[™]

Assured Gradient, Non-Elastic, Adjustable Therapy for Maintenance and Reduction of Lymphedema.

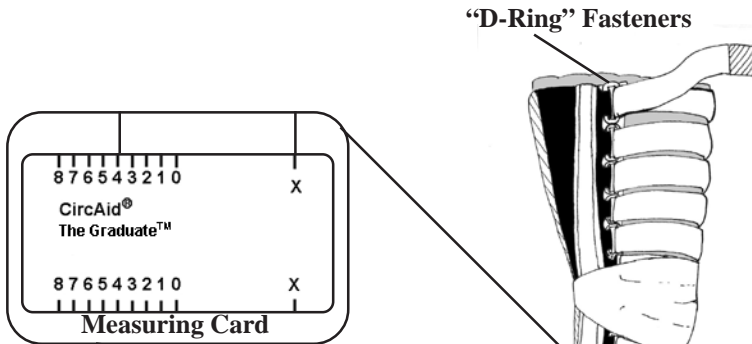

- ◆ Assured, Gradient Compression
- ◆ Non-elastic bands
- ◆ Immediately Adjustable
- ◆ Always Comfortable
- ◆ Gradient Counter Pressure
- ◆ Sustains Compression
- ◆ Flexible Knee
- ◆ Non-Elastic
- ◆ Night Time Garment
- ◆ Foam Liner Including Foot
- ◆ Minimal Adjustments
- ◆ One Year Warranty
- ◆ Hand Wash, Air Dry

Now you can be assured of a proper gradient with our new quick and easy Graduate[™] Thigh-High Legging. No more guessing. The Graduate features a series of non-elastic bands that tighten around a soft, comfortable foam liner. Through the use of "D-rings" and measuring guides, tightening and adjusting to gradient compression levels are easy to achieve. The foam liner insures smooth comfortable compression over the entire leg, **including the foot.**

An adjustable and flexible neoprene knee portion allows for an even greater range of leg motion. The Graduate is our most comfortable thigh-high legging yet.

The Graduate comes with an integrated foot piece which covers to the base of the toes. A strip of non-skid neoprene on the bottom of the foot provides an accommodating comfortable fit. This augments the legging in reducing and containing edema. The Graduate is available in custom sizes only due to various foot sizes.



A physicians prescription is required.