




Whole Leg EZ™

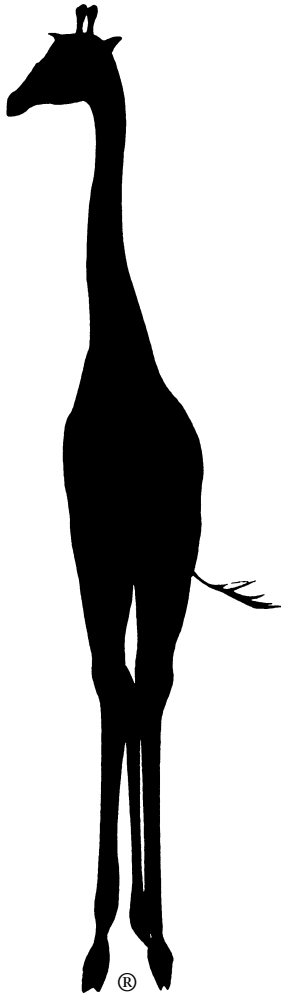


Our CircAid Systems have become the definition of clinically proven, non-elastic, adjustable and thereby sustainable compression in the treatment of venous and lymphatic disease.

- ◆ Easy On/Off
- ◆ Extremely Comfortable
- ◆ Flexible Knee
- ◆ Non-Elastic
- ◆ Immediately Adjustable
- ◆ Night Time Garment
- ◆ Adjustable Sustained Compression
- ◆ Foam Liner
- ◆ Gradient Counterpressure
- ◆ Minimal Adjustments
- ◆ Six Month Warranty
- ◆ Hand Wash, Air Dry

The CircAid Whole Leg EZ™ features a series of interlocking non-elastic bands that are tightened around a soft, comfortable foam liner. The foam liner insures smooth comfortable compression over the entire leg while also preventing the garment from sliding down the leg. A flexible knee portion of the garment allows for a greater range of motion when bending the leg at the knee. A strip of neoprene running along the back helps the garment conform to the shape of the limb. The Whole Leg EZ is easy to use and comfortable to wear.

The Whole Leg EZ comes with a CircPlus I Ankle-Foot Wrap which covers the ankle and foot area completely. This augments the legging in reducing and containing edema. The CircAid Thigh-High 2000 is available in standard sizes (small, medium, large, and X-large) which require eight circumferential and two length measurements. **Custom Whole Leg EZ leggings are also available.**

A physician's prescription is required.

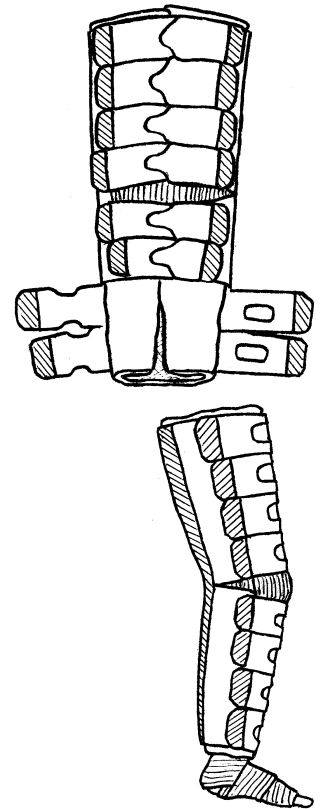

(rev0109)

Go to our website:
www.circaid.com
for full color pictures and short videos
Phone: (800) CIRCAID (247-2243)